


Freedom Wing Joystick Adapter

USER QUICK GUIDE

Introduction

The FreedomWing Joystick Adapter provides a way to connect a powerchair joystick to a gaming system. This can enable someone to game with the same joystick they use to control their powerchair. This adapter is compatible with joysticks that attach to the powerchair via a DB9 connector.

Features

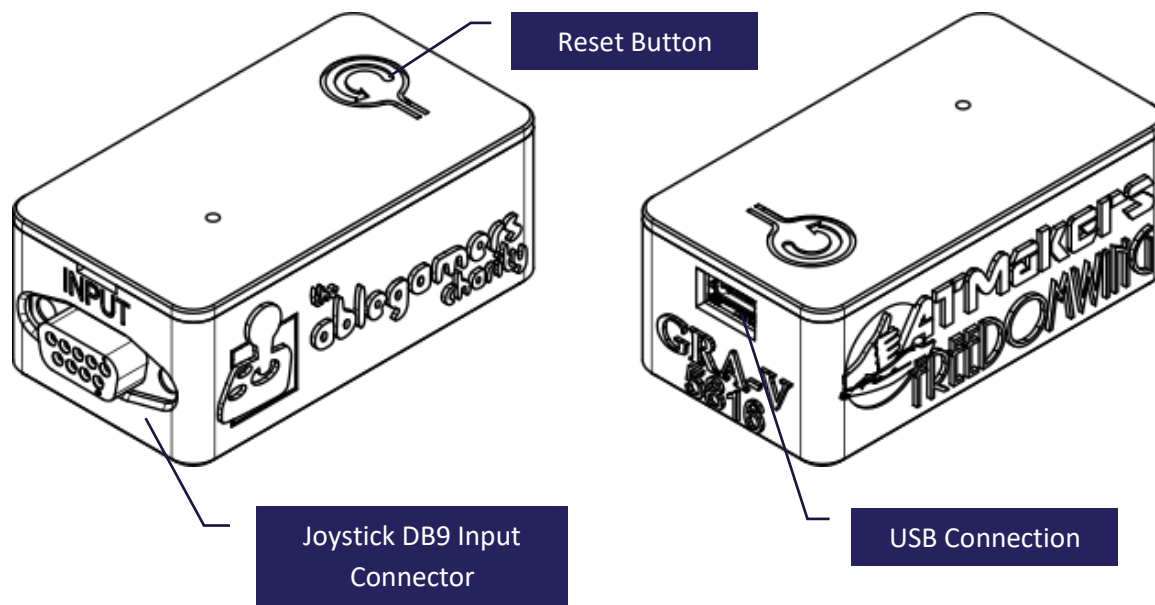

Setup Instructions

1. Disable wheelchair motors.
2. Disconnect the joystick DB9 connector from wheelchair.
3. Connect the joystick DB9 connector to the FreedomWing Adapter Input.
4. Connect the FreedomWing Adapter USB cable to the host device. (Use a suitable adapter if necessary.)

For use with the Xbox Adaptive Controller, it is recommended that the XAC be connected to a power source.

Freedom Wing Joystick Adapter

USER QUICK GUIDE

Compatibility

Wheelchair Manufacturer	Wheelchair Models	Required Interface
Merits Health USA	<ul style="list-style-type: none"> Merits vector P323 	<ul style="list-style-type: none"> R-Net OMNI interface R-Net OMNI2 interface
Permobil	<ul style="list-style-type: none"> F5 Corpus VS F5 Corpus F3 Corpus M5 Corpus M3 Corpus M300 Corpus HD M300 PS Jr. K300 PS Jr. Koola R-Net K450 MX F5 Corpus VS 	<ul style="list-style-type: none"> R-Net OMNI interface R-Net OMNI2 interface
Quickie	<ul style="list-style-type: none"> QM-710 QM-720 S636 Q400 M Q500 M 	<ul style="list-style-type: none"> R-Net OMNI interface R-Net OMNI2 interface

Cleaning

The FreedomWing Adapter can be wiped with a dampened cloth. Do not submerge the device.

Warnings

The powerchair should be completely inoperable while you are gaming.